


Actto 50



Actto 50

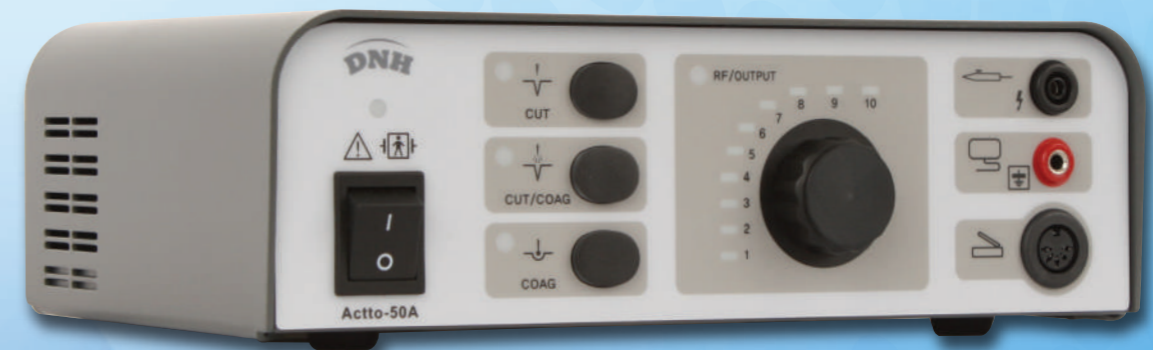

Included Components	QTY
Actto 50A Unit 220/230V	1
Power Cord	1
Reusable Electrode Set(6 type)	1
Dispersive Electrode with Cord	1
Footswitch assembly(Waterproof IP68)	1
Handpiece	1
Handpiece Stand	1
Owner's Guide	1

HIGH FREQUENCY THE COMPACT ELECTROSURGICAL UNIT

- + Exceptional cutting precision
- + Stop bleeding quickly
- + Fast healing

DNHmedtech Co.

Address : #211/F-dong, 70 Jimok-ro, Paju-si, Gyeonggi-do, Korea(10880)

Tel : (+82)70-4242-1011 Fax : (+82)2-318-2727


e-Mail : DNHmedtech@naver.com URL : www.DNHmedtech.com

CE
2460

DNH

DNH MedTech provides technology with improved performance at low power settings minimizing the risk and damage, and stop bleeding quickly

Actto 50

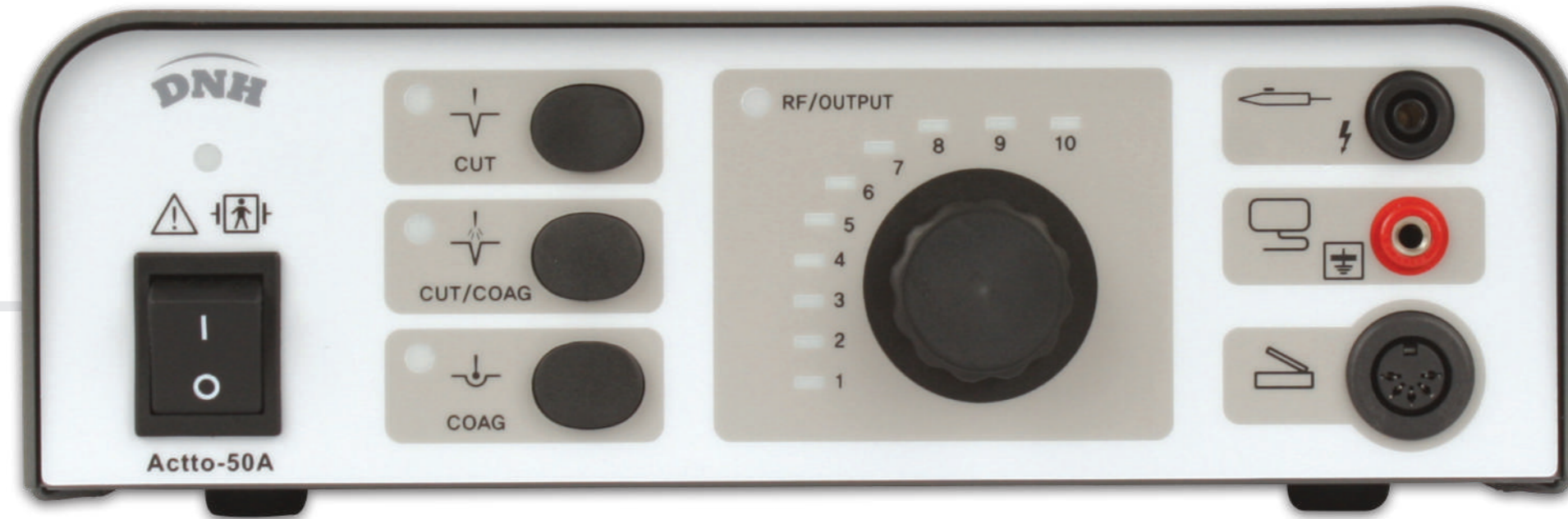

Actto 50 unit specifications

Width	10.4inches (265mm)
Height	3.5inches (89mm)
Depth	7.1inches (180mm)
Weight	6.8pounds (3.1kg)

Output Power (10 Levels)

CUT	50Watts
CUT/COAG	40Watts
COAG	30Watts

Power Supply 220/230V, 50/60Hz, 220VA

Frequency

1.2MHz ±10%

Electrical Classification class1, BF type

Actto 50, a high frequency electrosurgery system helps manage soft tissues with greater efficiency, greater precision and less bleeding than the common dental scalpel

- Remove unwanted tissue with ease
- Reduce chair time by simplifying operative procedures and maintaining hemostasis
- Provides outstanding cutting precision for superior clinical and aesthetic results
- Allows uneventful healing

Other clinical applications

- Widening gingival sulcus
- Crown lengthening
- Exposing impacted teeth
- Removing hyperplastic gingiva
- Frenectomy
- Stop bleeding quickly

Clinical Cases


• Gingivectomy

For extremely precise tracing of the ideal gingival line



• Exposure of impacted tooth

Circular incision for a clean and successful exposure

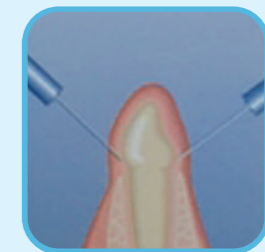


• Abscess incision and drainage

Fast and deep treatment



• Dental neck exposure

Ideal for treating hidden decay

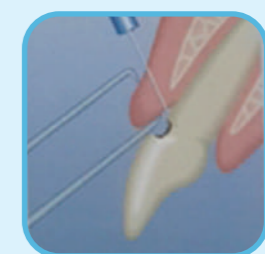

Actto 50